

# WA-900 SUPERVISED 900MHz Wireless WANDERER ALERT™

**WORKS AUTOMATICALLY!**

**Detects patients wandering away and then sounds the alarm and calls for help!**




**1 Transmitter Only Per  
WA-900 System.**

**NOW YOU CAN PROTECT YOUR LOVED ONES AND  
PATIENTS FROM WANDERING AWAY!**



## **WANDERER ALERT™**

The pendant Transmitter automatically sends a signal every 10 seconds to the receiver-dialer. If the receiver doesn't hear from the pendant within an adjustable period, (indicating patient has traveled a distance away) the receiver-dialer will automatically sound the alarm strobe and/or dial emergency numbers and deliver the emergency message.

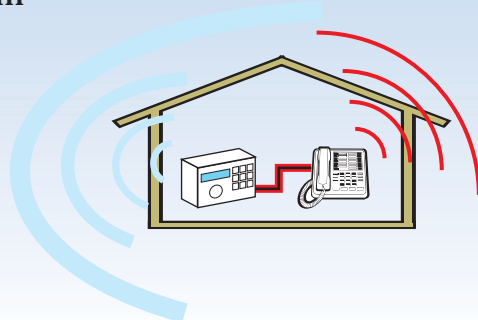

## WHEN ACTIVATED...

The Wanderer Alert automatically detects wandering, sounds the alarm and calls to deliver the message for help to:



+ Family or neighbors



+ Emergency services  
(local regulations permitting)




+ Health Care Providers



+ Cellular phones and pagers

## Wireless 900 MHz Technology comes to the rescue of loved ones and patients

### How It Works:

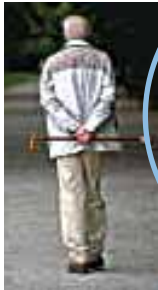

**1.** The wireless pendant is attached to the patient with a non-removable band and sends a continuous signal to the Wanderer Alert unit every 10 seconds.



**2.** When the patient attempts to wander away, the Wanderer Alert unit detects the missing signals and sounds the horn/strobe alarm and/or dials the telephone, cell or pager numbers that you have installed and leaves your message in any language. A switch on unit allows you to select dial out.

**3.** The Wanderer Alert system is capable of monitoring both High Risk and Low Risk patients.

Wanderer Alert Automatically Sounds the Horn Strobe and calls for help, even if care provider is away! **No Monitoring Fees!**

### High Risk Patients

Patients that are high risk to wander away.

For High Risk patients, turn Range Dial Left or CCW until desired Safe Area or Roaming Range is satisfactory. The pendant will activate Wanderer Alert when patient wanders from unit, allowing patient to move throughout home or living area.

### Low Risk Patients

Patients that are low risk to wander away.

With Low Risk patient, turn Range Dial Right or CW until desired Safe Area or Roaming Range is satisfactory. This allows them to move throughout living areas and into garage, and areas such as the yard, garden and to the mailbox. Wanderer Alert will activate when patient wanders from the unit.

Pendant may be Worn on Wrist or Ankle, Using Supplied Band.



**1 Transmitter Only Per WA-900 System.**

### Order Information

| Model          | Description   |
|----------------|---|
| WA-900 .....   | Complete with Wanderer Alert Unit, Horn-Strobe, wireless Pendant with extra wrist ankle locking bands, AC Adapter |
| PP-1 .....     | 24 Hour Battery Back Up Pack  |
| EM-900-WA .... | Monitors up to 40 different patients  |
| LWA-900 .....  | Patient Locator   |

See Website Product Index and Price List for all Accessories for Products and Systems.

## F E A T U R E S

- + Sends a prerecorded emergency voice message (up to 20 seconds in length) three times to four different numbers. Allows for "listen in" and "2-Way" speakerphone communication mode.
- + No installation required! Just plug into any standard telephone jack, program the emergency telephone numbers to dial and record voice message in any language requesting help.
- + Automatically activated by pendant transmitter when patient wanders away.
- + Optional 24-hour battery back-up ensures operation even when the power fails.
- + The WA-900™ automatically skips to the next number when a busy signal is detected or when telephone is not answered. The WA-900™ then keeps calling back until message is delivered or The WA-900™ is deactivated.



Contact Us: (T) 800-227-2474 (F) 800-2271592

www.unitedsecurity.com  
www.emergencycaller.com